

1851 East Paradise Road, Suite A
Tracy, CA 95304
Phone: 866.592.0488 • Fax: 209.834.8832
topshelforthopedics.com



Solar Walker

Instructions for Use

Please read these instructions completely and carefully before applying this product. Correct application is important for proper function of this brace.




Solar Walker

Instructions for Use



Please read these instructions completely and carefully before applying this product. Correct application is important for proper function of this brace.

Indications: Soft tissue injuries of the lower leg, acute ankle sprains, stress fractures of the lower leg, stable fractures of the foot and ankle, stable delayed union or non-union fractures of distal tibia and fibula, achilles tendon repairs.



Application Instructions:

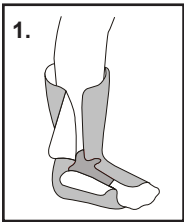

1. For first time application - loosen straps and remove liner from the boot.
2. Place foot in the liner with the heel fitting snugly in the posterior section and secure the foot portion then the calf portion with the contact closure.
3. While spreading the metal uprights on the boot, position the foot with the liner affixed into the base of the boot aligning the uprights to the midline of the ankle. Remove the plastic sleeves covering the metal uprights to secure the liner to the walker struts.
4. Secure the hook/loop straps on the front of the boot walker starting at the toes and working up the brace towards the calf.

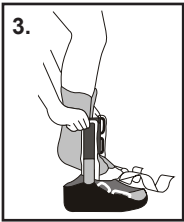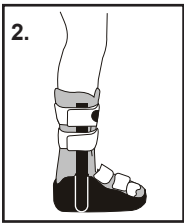

Internal Air Pressure Control: To put air inside the boot, turn the valve clockwise until secure, then squeeze the bulb until you feel adequate pressure. To remove air from the boot, simply turn the valve counter clockwise.

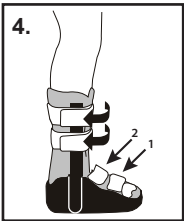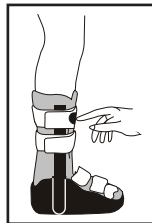

Removal: To remove the walker boot open completely the front straps on the boot walker and open the liner closure. Do not remove the liner from the struts.

Re-Application: When re-applying the walker after the initial fitting it is not necessary to remove the liner from the boot, simply re-secure the liner closure and lace the straps from the foot to the calf.

Use and Care:

Hand wash the liner with mild soap and water and lay out to dry. Do NOT machine dry.

Notes:

For single patient use only.

Warning:

This product is to be used under the supervision of a medical professional. If you experience any pain, swelling, sensation changes, or any unusual reactions while using this product, consult your medical professional immediately.

Warranty:

Top Shelf Orthopedics guarantees this product to be free of defects in material and construction for a period of six months from the date of purchase.

Caution:

Federal Law (U.S.A.) restricts this device for sale by or on the order of a licensed healthcare professional.